

To ORDER

Ref.	Use	Length
NBI035	With rigid cystoscope	35 cm
NBI070	With flexible cystoscope	70 cm

INDIVIDUALLY PACKED - STERILE - 

BoNee

Needle for bladder wall injection




BoNee

Locking

The Touhy-Borst connector allows to lock the needle to the cystoscope working channel

- Facilitates the puncture
- Facilitates needle rotation

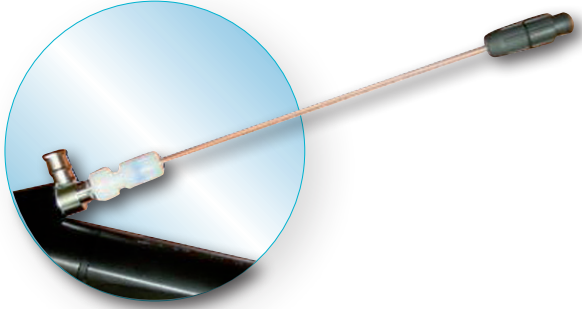


Designed with urologists to facilitate bladder wall injection

Can be locked →



6 Ch/Fr flexible and strong body

- Compatible with flexible and rigid cystoscopes
- Allows to perform in any part of the bladder without kinking

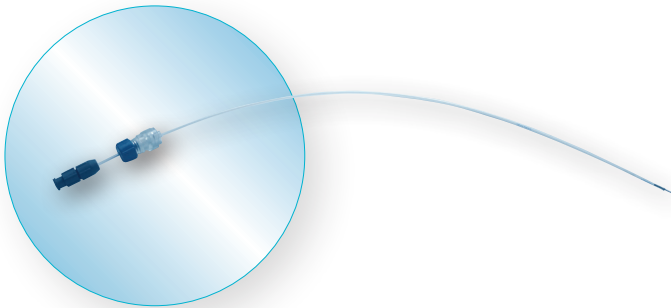

6 Ch/Fr (1,7mm) flexible and strong body →

4mm tip

WITH PROTECTION CAP

- Injection inside the bladder muscle
- The protection cap avoids damaging the working channel of the cystoscope during needle insertion*



* Please remove the protection cap before inserting the cystoscope in the urethra - see IFUs

4 mm tip - Chiba 22 G with protection cap →

