


To ORDER

Ref.	Use	Box	Length
BBI035	With rigid cystoscope	5 pcs	35 cm
BBI070	With flexible cystoscope	5 pcs	70 cm

INDIVIDUALLY PACKED - STERILE - 

BoNee

Needle for bladder wall injection

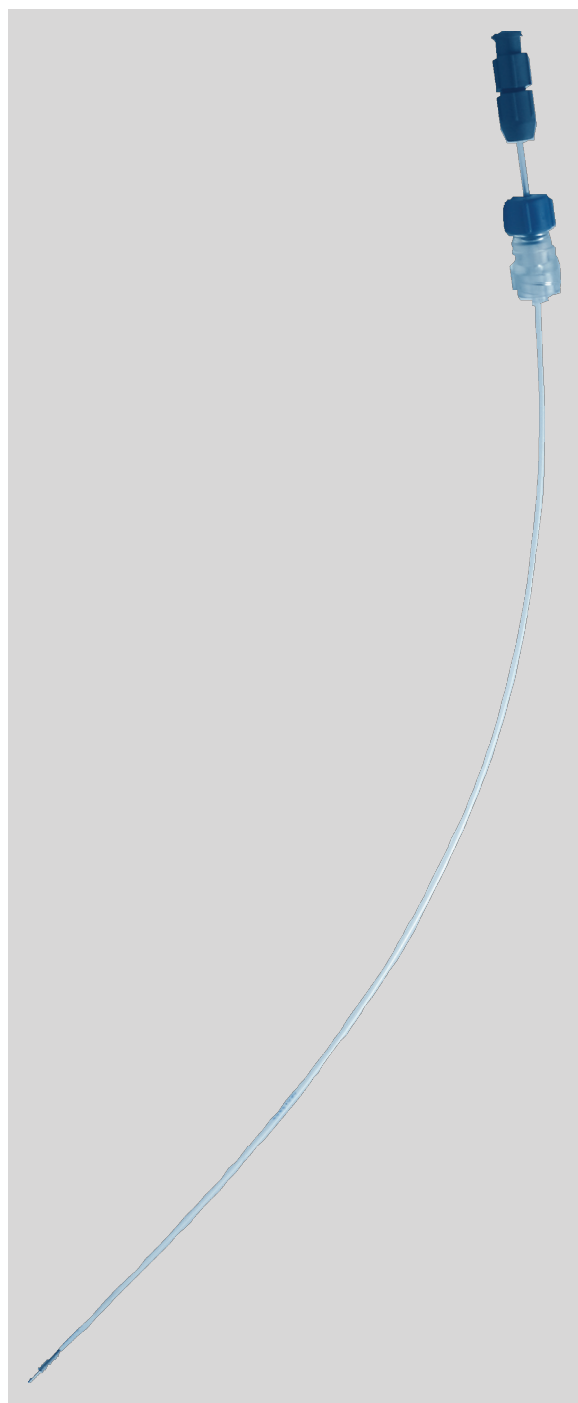


BoNee

Locking


The Touhy-Borst connector allows to lock the needle to the cystoscope working channel

- Facilitates the puncture
- Facilitates needle rotation

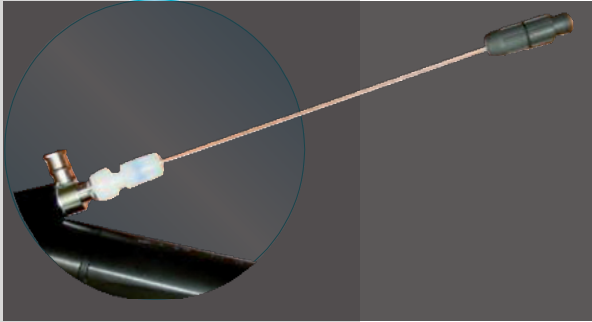

6 Ch/Fr flexible and strong body


- Compatible with flexible and rigid cystoscopes
- Allows to perform in any part of the bladder without kinking

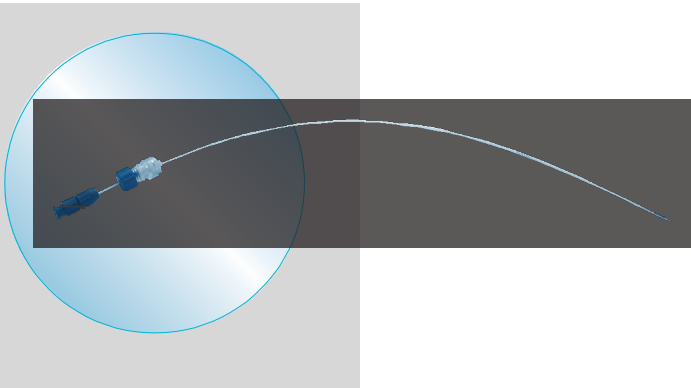

4mm tip

WITH PROTECTION CAP

- Injection inside the bladder muscle
- The protection cap avoids damaging the working channel of the cystoscope during needle insertion*

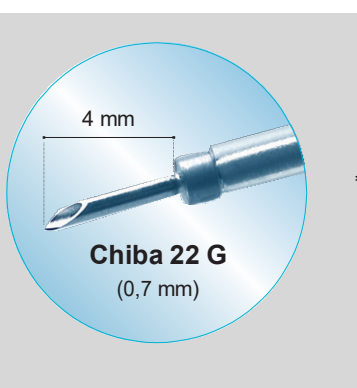


* Please remove the protection cap before inserting the cystoscope in the urethra - see IFUs

Designed with urologists to facilitate bladder wall injection

Can be locked →

6 Ch/Fr (1,7mm) flexible and strong body →

4 mm tip - Chiba 22 G with protection cap →

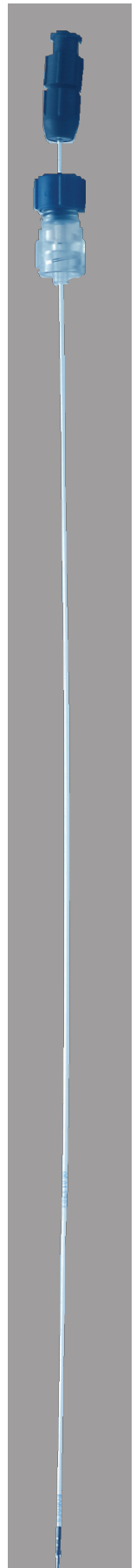

* Developed in partnership with Dr. Corcos, McGill University, Montreal, Canada